

Duplicating Dentures

THIS TECHNIQUE WILL ACCURATELY RECORD THREE DIMENSIONAL INFORMATION

- **RECORDS VERTICAL DIMENSION**
- **ESTABLISHES BITE REGISTRATION**
- **AIDS IN TOOTH SELECTION**

OTHER BENEFITS INCLUDE:

- **REDUCED CHAIR TIME**
- **VERY SIMPLE AND FAST**
- **PATIENT DOES NOT HAVE TO BE
WITHOUT DENTURE**

Lab Putty



◆ **EXAFLEX PUTTY**
BY
GC AMERICA

◆ **800-323-7063**



◆ **SIL TECH SUPER**
BY
IVOCLAR

◆ **800-533-6825**

(Sil-Tech is more rigid, but both will work.)

Putty Method




- ◆ **BEGIN BY MIXING 4 SPOONS OF A WITH 4 SPOONS OF B.**

Putty Method



- ◆ PRESS AND ADAPT INSIDE TISSUE AREA OF DENTURE.
- ◆ FORM A BASE ON THE COUNTER & SET DENTURE ON BASE.
- ◆ CHECK PERIPHERIES FOR GOOD ADAPTATION.
- ◆ LET SET UP ABOUT 5 MINUTES. **DO NOT SEPARATE**

Putty Method

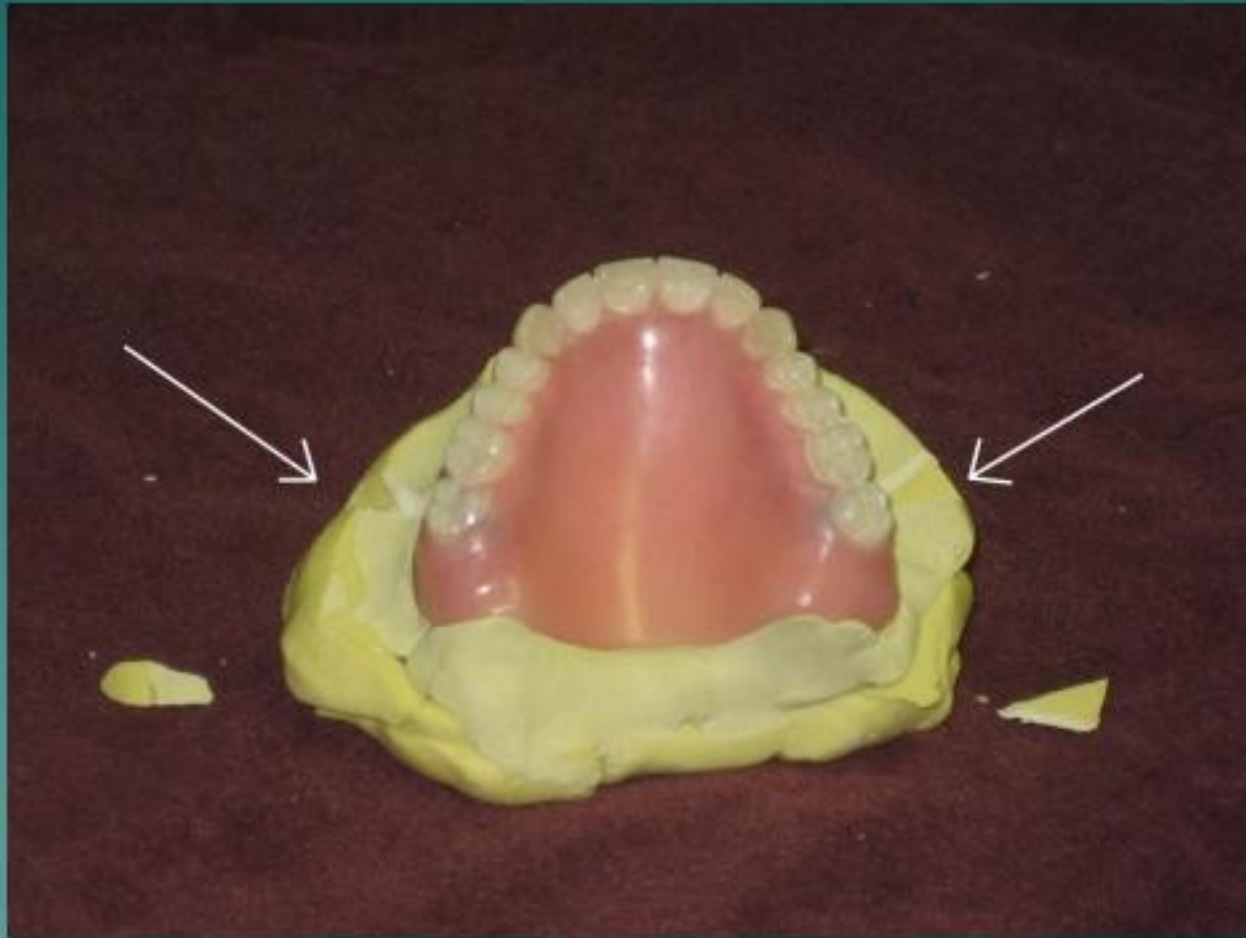

- ◆ **CUT KEY NOTCHES IN OUTSIDE BASE.**

Putty Method



- ◆ LIGHTLY VASELINE OUTSIDE OF BASE SURFACE.

Putty Method



- ◆ **MIX 4 SPOONS OF A AND 4 SPOONS OF B.**
- ◆ **ADAPT TO TOP OF DENTURE TEETH AND SURFACE.**

Putty Method



- ◆ **ADD TO FORM A MOUND AGAINST THE BASE.**
- ◆ **ALLOW 5 MINUTES TO SET UP.**

Putty Method



- ◆ **SEPARATE AND INSPECT DUPLICATION.**
- ◆ **DUPLICATE OTHER DENTURE *(IF APPLICABLE)*.**
- ◆ **RECORD SHADE AND RETURN DENTURE TO PATIENT.**

Putty Method

- ◆ **SEND TO LABORATORY WITH THE FINAL IMPRESSION, OR A WASH CAN BE DONE IN THE RETURNED SET-UPS AT NEXT APPOINTMENT.**